

PERFECT FIT

COLLECTION

STERILIZATION CENTER

Hello perfect work flow.



MCC
Modular & Custom Cabinets



TAILORED SOLUTIONS

FOR YOUR DENTAL PRACTICE

Our Perfect Fit Sterilization centers are designed to support the most efficient work flow for your practice's sterilization needs.



'Green' construction. Our Perfect Fit Cabinetry is also available using environmentally friendly manufacturing techniques. 'Green' construction means our cabinets contain no urea formaldehyde resins and are built with water-based adhesives. As a result, they are LEEDS compliant.

This new design is the result of our collaborative process, incorporating ideas from dental practitioners across the country.

- The collective goal was to streamline the flow of the sterilization process, aid operator efficiency and surpass the latest CDC asepsis guidelines.
- But at MCC, we've taken our commitment even further - with a strong emphasis on 'green' materials and construction.
- Wherever possible we utilize solid, non-porous surfaces, that emit no urea formaldehyde fumes and absorb no moisture.
- The result is a unit that satisfies all the criteria necessary to make it LEEDS compliant. This 'Leadership in Energy and Environmental Design' certification means the sterilization center is a safe choice for your staff, your patients and the environment, for years to come.

Investing in well designed, high-quality cabinetry today pays high dividends down the road, as practices evolve to meet new and ever-demanding needs.

Talk to our experts about any special features you may require. A few simple custom modifications are a cost-effective way to guarantee the 'perfect fit' for your practice.



RECEIVING & CLEANING



PREPARATION & PACKAGING



STERILIZATION



STERILE CLEAN STORAGE

STAGE 1

RECEIVING & CLEANING

- 5-1/2 feet of feature-packed space has been allocated to this stage.
- Ample room for a cassette-based environment
- Use of cassettes requires 60% less linear counter space making it perfect for compact installations
- Upper Red Light Module is designed for dirty instrument storage until processing is initiated
- Shelves are constructed from sturdy & moisture-resistant 1/4" solid surface material
- Storage area for a hand piece maintenance system has been incorporated into the layout
- Solid surface wrap drawer covers protect drawer contents and provide the operator with additional work surface
- Standard hands-free water activation, large sink and goose neck faucet accommodate oversized items
- Integrated waste and sharps disposal
- Convenient power access panel for additional electrical requirements
- Use of tabletop/built-in ultrasonic or instrument washer



STAGE 2

PREPARATION & PACKAGING

- Side barriers protect materials during this key stage of the process
- Barriers are made from a high density/ non-porous material that will not rot, splinter, swell or delaminate
- Our ST Sterilization Centers are built with a high-strength scratch resistant surface which inhibits growth of damaging microorganisms.
- Barriers combined with the solid surface countertop form the perfect anti-microbial area in which to package and process materials
- Solid surface wrap drawer covers protect drawer contents and provide the operator with additional work surface
- Lower cabinets with pull-outs assist personnel in accessing sterilization supplies

STAGE
3

STERILIZATION

- Sterilizers have been located with ergonomics in mind
- Series of three slide-out counters in the tower section accommodate most sterilizers and Statim while providing a central slide tray/work surface for convenient loading and unloading
- Lower drawer is ideal for the storage of sterilization materials
- Construction of this tower area includes similar high density/non-porous materials utilized throughout
- An exhaust fan has been installed above the sterilizer in the unique rear chimney area of the tower, providing the proper operating temperature and preventing expensive equipment repairs
- To assist in installation, an access panel has been built into the chimney area



STAGE
4

STERILE CLEAN STORAGE

- Upper storage module provides sufficient storage for trays and cassettes
- Upper storage includes Blue Light Module to signify items have been processed and are ready for use
- Shelves are constructed from sturdy, moisture-resistant 1/4" solid surface material routed to resist movement on shelf supports
- An ample number of storage drawers are provided below the counter with one drawer above



FEATURED COLORS

Choose from a wide variety of the latest solid colors and designer hues and classic wood grain patterns.

LAMINATE STANDARD COLLECTION



White • 949-58



Antique White • 932-58



Folkstone • 927-58



Crystal • D388-60 (Surf)



Prestige Maple • L476 Matte Finish
(Amber Curly Maple)



Pearwood • WX031



Versailles Anigre • 7923-60
(Slipping Seattle Java)




Biltmore Cherry • 7924-60



Tuxedo • L477



Starry Night • L762 Matte Finish

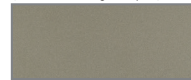

Milky Way • L767 Matte Finish



Portobello • L495 Urbania Finish



Truffle • L479 Matte Finish



Summer Breeze • L540 Alto Finish



White Chocolate • L492 Urbania Finish



Latitude North • L530 Alto Finish




Maritime Dune • L555 Urbania Finish



Brushed Aluminum • L474 Matte Finish

QUARTZ COLLECTION by caesarstone Quartz Surfaces

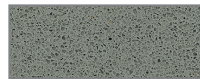

Cement • 3040



Jet Black • 3100




Stone Grey • 4030




Creme Brule • 4255



Ocean Foam • 6141

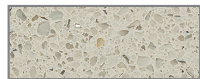

Quartz Reflections • 7141

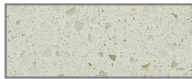


Ice Snow • 9141

SOLID SURFACE COLLECTION




Light Beige • 254

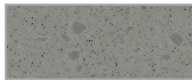


Blue Bayou • 510




Arctic Granite • 603

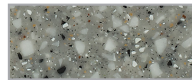


Platinum Granite • 615

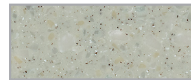

Moonstone Granite • 674

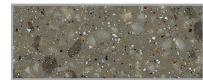


Cinder Granite • 679 AR

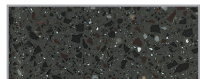


Midnight Granite • 690



Sanibel Granite • 691

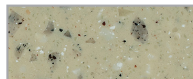

Rocky Road Granite • 729 AR

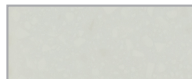


Aspen Granite • 780

Let us help you make your color choices.

We're here to help you choose the surface colors that best match your office style and design objectives. Ask us for assistance. We'll be happy to show you all the color choices available to you.

Other colors are available upon request, up-charge will apply. Prices and colors are subject to change without notice. The products and colors shown in this brochure should only be used as a guide and are not absolute and may vary. Please contact one of our sales representatives for accurate color samples and information.

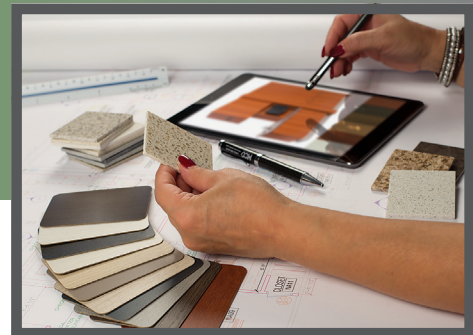

**TAILORED SOLUTIONS
FOR YOUR DENTAL PRACTICE**

Sterilization Centers • Rear Cabinets • Side Cabinets • Center Island Units • Mobiles • Orthodontic Mobiles • Pediatric Benches • Brilliance Units

1-800-388-6236

www.mccdental.com

FOLLOW US ON THE MAJOR SOCIAL NETWORKS

